

HISTOACRYL® CLASSIC

HISTOACRYL® FLEXIBLE



6501103V
Histoacryl® Clear



6501102V
Histoacryl® Blue

Applicator Tip



32028864
Histoacryl® Flexible

Histoacryl® Flexible directly compares to Dermabond Advanced

Histoacryl® Tissue Seal Comparison

	Histoacryl® Classic	Histoacryl® Flexible
Composition	n-butyl-cyanoacrylate	n-butyl-cyanoacrylate Plasticizers
Color	<ul style="list-style-type: none"> • Violet (D&C violet n°2) • Clear 	<ul style="list-style-type: none"> • Violet (D&C violet n°2)
Indications	Skin Closure	Skin Closure
Application	Direct from ampoule	Direct from ampoule
Color Code (ampoule)	Translucent & Yellow bar	Opaque ampoule & Yellow bar
Color Code (packaging)	White and Orange	Dark and Light Blue
Storage	24 months @ 22°C	24 months @ 25°C

Why Histoacryl® Flexible?

- Histoacryl® Flexible has been designed to combine flexibility and strength to close and protect wounds
- Peel strength is significantly higher (p. O.OO1) for Histoacryl® Flexible than for Dermabond
- Flexibility allows closure of incisions up to 25cm
- Only one layer of adhesive is necessary to achieve a full-strength bond
- Special applicator tip provides better control and improved application
- Cyanoacrylate-based topical skin adhesives yield excellent cosmetic results
- Generates less heat during polymerization, making it more comfortable for the patient
- Adapts easier to the contours of the surgical incision

HISTOACRYL® TISSUE SEAL

Topical Skin Adhesive

Microbial Barrier

In vitro studies have shown that Histoacryl® Flexible acts as a barrier to microbial penetration as long as the adhesive film remains intact. Histoacryl® Flexible inhibits the growth of Staphylococcus aureus, Staphylococcus epidermidis, Escherichia coli, Pseudomonas aeruginosa, Enterococcus faecium, Brevundimonas dimuta and Candida albicans for up to 7 days on inoculated agar plates.

Part #	Cross Ref.	Description	UOM	#/UOM
32028864	TS1051260P	Histoacryl® Flexible 10 / 0.5 mL ampoules 10 Applicator Tips	Box	10
6501102V	TS1050044FP	Histoacryl® Classic - Blue 10 / 0.5 mL ampoules 10 Applicator Tips	Box	10
6501103V	TS1050071FP	Histoacryl® Classic - Clear 10 / 0.5 mL ampoules 10 Applicator Tips	Box	10

OPEN



Open the package and take out the Histoacryl® Tip.

PREPARE APPLICATOR TIP



Open the ampoule by twisting off the ribbed tip. Attach the Histoacryl® Tip to opened ampoule.

APPLY




Apply the glue in a thin layer to the approximated wound edges and hold in apposition for 60 seconds.

The polymerization starts immediately upon contact with tissue.

PACKAGING



Box
10 Ampoules
10 Applicator Tips

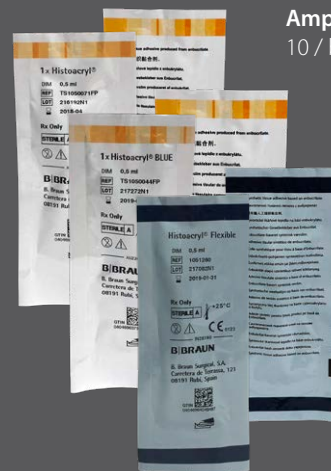

Ampoules
10 / box



Applicator Tips
10 / box